


SilFlow™ Deans' Switch Kit

Installation Manual

1. INTRODUCTION

The SilFlow™ Deans' Switch is a new chromatography solution using SGE's proprietary SilFlow™ microchannel technology.

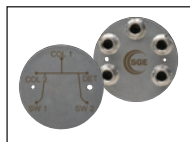

SilFlow™ Deans' Switch

SilFlow™ is an innovation in design and fabrication resulting in a highly efficient and reliable micro fluidic platform. The SilFlow™ Deans' Switch is based on the famous Dean's principle, incorporating advanced technical improvements. The Deans' Switch forms the core of most multidimensional chromatography configurations.

SilFlow™ Deans' Switch Kit is manufactured as a capillary column gas flow switching system. Applications are numerous such as heart cutting, column back flushing, inlet switching, column switching, detector switching, solvent venting etc.

Features of SilFlow™ Deans' Switch

- Chemically deactivated internal channels, resulting in analysis with enhanced quantitative accuracy and high reproducibility. Can be incorporated into your system without impacting your chromatography.
- Low dead volume connections, SilTite™ FingerTite metal ferrules result in a reliable zero dead volume connection, giving you optimized peak shapes.
- Superb operational stability. Thermal Lag – tracks oven temperature up to 20 °C/min, the design alleviates cold spots and sample condensation. Maximum temperature – No practical temperature limit. Limited only by the temperature rating of the GC column being used (420 °C).
- Easy to install and leak free. The kit incorporates SilTite™ FingerTite fittings that are easy to set up and can be tightened

using finger force to achieve a reliable seal, even for the most sensitive MS systems – no wrenches are required.

2. CONTENTS OF THE SILFLOW™ DEANS' SWITCH KIT

SilFlow™ Deans' Switch Kit contains accessory parts, tools and the installation manual. Parts and tools in the Kit are listed in Table 1 and shown in Figure 1.

Part Number	Part Description
123726	SilFlow™, Deans' Switch
-	Mounting Bracket
123751	2 m x 100 µm/363 µm VSD Tubing
123755	2 x 1.1 mm OD SS Tubing with 1/16" sleeve X 1 m
123706	SilFlow™ SilTite™ Ferrule 0.4 mm, Packet of 10
123707	SilFlow™ SilTite™ Ferrule 0.5 mm, Packet of 10
123705	SilFlow™ SilTite™ Ferrule 1.1 mm, Packet of 10
123704	SilTite™ FingerTite Nuts, Packet of 6
-	Screw for Mounting Bracket
123717	FingerTite Tool
-	Pre Swage Tool, 0.4 mm
-	Pre Swage Tool, 0.5 mm
-	Blanking Pins, Pack of 2

Table 1. Parts and Tools in the SilFlow™ Deans' Switch Kit.

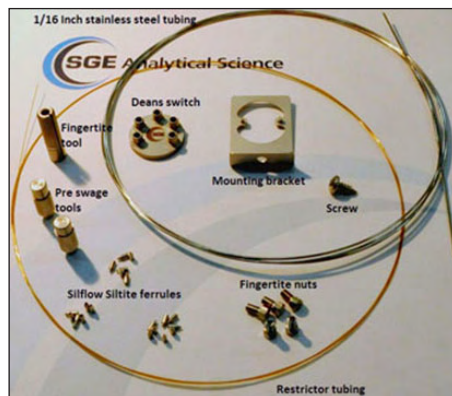


Figure 1. Kit contents – Parts and Tools identification.

3. INSTALLING THE SILFLOW™ DEANS' SWITCH KIT

In these instructions it has been assumed the configuration is for heart cutting two dimensional gas chromatography applications (Figure 2). Please email support@sge.com for specific application requirements.

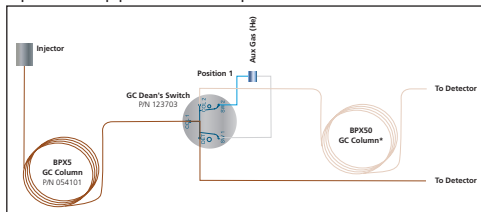


Figure 2. SilFlow™ Deans' Switch set-up for heart cutting two dimensional gas chromatography.

3.1 Installing the Mounting Bracket

Prior to installation ensure the GC oven is turned off and cooled. Installing the mounting bracket is straight forward. Make sure the interior of the oven is free from obstacles such as GC columns and column holders as this will make the installation easier.

- i. Identify a convenient location in the GC oven. Consider the fact that after installation it should be easy to access the components of the system: columns, detectors, injectors and accessories.
- ii. Find an appropriate hole in the location of the oven wall or carefully drill a 3.2 mm sized hole in order to attach the screw.
- iii. Secure the bracket against oven wall and align the screw to the fitting and into the wall (Figure 3).

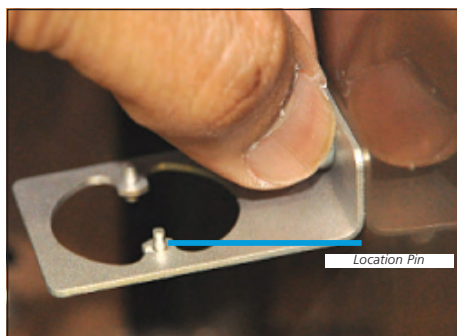


Figure 3. Attaching the mounting bracket.

- iv. Tighten the screw carefully with a screwdriver.

3.2 Attaching the SilFlow™ Deans Switch

- i. Mount the SilFlow™ Deans' Switch to the oven bracket with the guidance of the locator pins ensuring it sits on the bracket firmly (Figure 4).

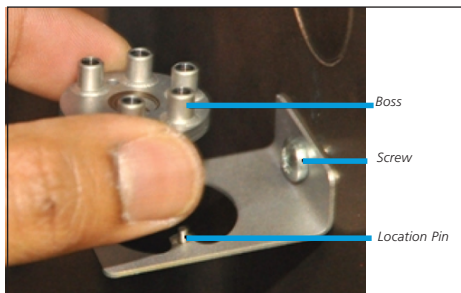


Figure 4. Attaching the Deans switch.

Make sure the bosses face the correct direction for easy access to install the tubing.

- ii. Identify the appropriate port locations marked on the back of the SilFlow™ Deans' Switch (Figure 5). Locate the appropriate SilFlow™ ferrules and nuts and prepare to connect the relevant tubing to the SilFlow™ Deans' Switch.

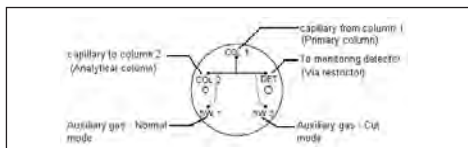


Figure 5. Identifying Deans switch port locations.

3.3 Installing the Columns

- i. Re-install any column holders previously removed from the oven in order to attach the mounting bracket.
- ii. Hook the GC columns (1st and 2nd Dimension) onto the column holders firmly.
- iii. Install one end of the 1st Dimension GC column (Primary) to the injector in accordance to the GC manufacturer's recommendation.
- iv. Connect the other end of the 1st Dimension GC column to the 'COL 1' port of the Deans switch.
- v. Likewise install the 2nd Dimension GC column between the 'COL 2' port of the SilFlow™ Deans' Switch and the detector system (2nd Dimension).

When connecting the capillary column to SilFlow™ Deans' switch, use the Silflow™ FingerTite Jig to pre-swage the ferrule (Figure 6) before attaching the column to the SilFlow™ Deans' Switch boss,

ensuring the capillary will not be crushed onto the SilFlow™ Deans' Switch channels.

There are two types of jigs available; one for the 0.4 mm ID ferrules (for 0.36 OD capillaries) and one for the 0.5 mm ID ferrules (for 0.43 OD capillaries).

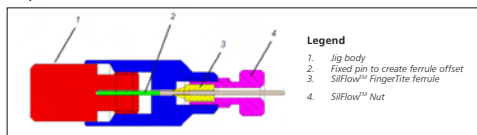


Figure 6. Assembly with SilFlow™ FingerTite Jig.

Use the FingerTite tool to tighten the FingerTite nuts (Figure 7) applying finger force only. Never use a wrench.

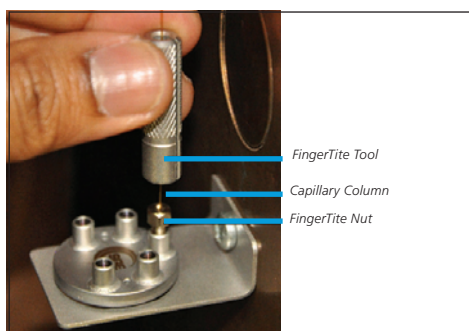


Figure 7. Using the FingerTite Tool.

3.4 Installing the Restrictor Tubing

The dimensions of the restrictor tubing should be carefully selected in such a way that the gas flow resistance is similar to that of the 2nd Dimension GC column (Analytical). To keep the time lag short as possible for the analytes to reach the 1st Dimension detector (Monitoring) the restrictor has to be short in length.

The balancing restrictor can be constructed from the narrow bore uncoated deactivated fused silica tubing (VSD tubing) supplied with the kit. Use 0.1 mm ID tubing for the GC system with 0.22 – 0.25 mm ID columns and 0.15, 0.22 or 0.25 mm ID tubing for the larger diameter columns.

i. Calculate the dimension of the restrictor tube based on the operating pressures of the GC system and dimensions of the 2nd Dimension GC column. An example is shown in Table 2, restrictor tube has to be accurately trimmed until the flow rate matches that through the 2nd Dimension GC column.

- ii. Connect one end of the restrictor to the 'DET' port of the Deans' switch. Use the FingerTite Jig to pre-swage the ferrule.
- iii. Install the other end of the restrictor to the 1st Dimension detector system (Monitoring) in accordance to the GC manufacturer's recommendation.

For operating pressures of the GC system, refer to the recommendations of the column manufacturer and consider method conditions.

As an example, method parameters for a heart cutting two dimensional gas chromatography technique using BPX5 and BPX70 columns for the analysis of essential oil are listed in Table 2.

1 st Dimension GC column	30 m in length, 0.25 mm ID, 0.25 µm film; BPX5
2 nd Dimension GC column	30 m in length, 0.25 mm ID, 0.25 µm film; BPX70
Restrictor	0.1 mm ID VSD tubing 0.75 m in length
Carrier	Helium
Inlet Pressure	41.7 Psi
Mid point (Auxiliary) Pressure	27.5 Psi
Inlet Temperature	250 °C
Oven Temperature	140 °C
Detectors	FID at 300 °C

Table 2. Method parameters for a heart cutting two dimensional GC chromatography – Analysis of essential oil.

3.5 Installing the Stainless Steel Tubing for Auxiliary gas

- i. Install 1.1 mm OD end of the 1/16 inch sleeved stainless steel tube to the 'SW 1' port of the SilFlow™ Deans' Switch. Use the SilFlow™ Ferrule 1.1 mm and a FingerTite nut for the connection.
- ii. Likewise connect one end of the other 1/16 inch sleeved stainless steel tube to the 'SW 2' port of the SilFlow™ Deans' Switch.
- iii. Run the stainless steel tubes through the oven ceiling and exit to the location where the solenoid valve is installed. If your GC is not equipped with EPC, SGE recommends VSO-EP™ Electronic Pressure Control Unit. See www.parker.com/precisionfluidics. Bend the tubes and align in such a way that they run closely against the oven walls. This provides easy access to the oven (Figure 8).
- iv. Connect 1.1 mm OD end of the stainless steel tubes to the two outputs of the solenoid valve using the appropriate connection. Refer to the manufacturer's user manual for correct installation.

4. TAKING CARE OF YOUR SILFLOW™ DEANS' SWITCH SYSTEM

Tubing to the SilFlow™ Deans' Switch can be disconnected and reconnected many times without removing the pre-swaged ferrules. However, it is very important to inspect the capillary end carefully before reconnecting to the SilFlow™ Deans' Switch. Make sure the end is not squashed or cracked. If the capillary does not have a clean square end, the column end has to be cut again and a new ferrule to be pre-swaged.

When the SilFlow™ Deans' Switch is not connected, protect it from particulates and dust contaminating the internal channels. Make sure to block the bosses using either self-sealing nuts or pre-swaged ferrule to pieces of metal wire and appropriate SilFlow™ FingerTite nuts.

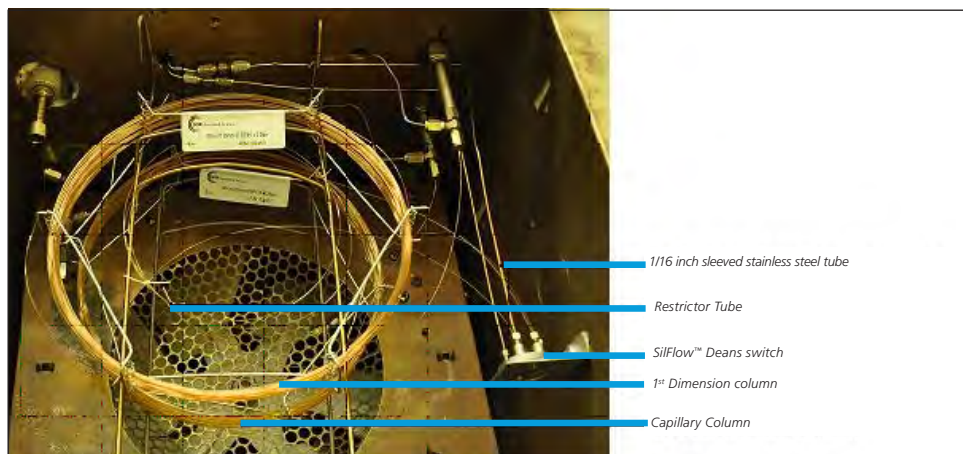


Figure 8. Hardware set-up for heart cutting two dimensional gas chromatography.

<p>AUSTRALIA & PACIFIC REGION SGE Analytical Science Pty Ltd Toll Free: 1800 800 167 Tel: +61 (0) 3 9837 4200 Fax: +61 (0) 3 9874 5672 Email: support@sge.com</p>	<p>CHINA SGE Shanghai Representative Office Tel: +86 21 6407 9382 Fax: +86 21 6407 9386 Email: china@sge.com</p>	<p>MIDDLE EAST SGE Gulf Tel: +971 6 557 3341 Fax: +971 6 557 3541 Email: gulfsupport@sge.com</p>	
<p>EUROPE SGE Europe Ltd European Head Office Toll Free: 00800 2790 8999 Toll Free Fax: 00800 2626 2609 Tel: +44 1908 568 844 Fax: +44 1908 566 790 Tel France: +33 1 69 29 80 90 Fax France: +33 1 69 29 09 25 Tel Germany: +49 (0) 6155 / 60746 0 Fax Germany: +49 (0) 6155 / 60746 0 Email: europe@sge.com</p>	<p>INDIA SGE Laboratory Accessories Pvt Ltd Tel: +91 22 24715896 Fax: +91 22 24716592 Email: sgeindia@vsnl.com</p>	<p>UNITED STATES OF AMERICA SGE Incorporated Toll Free: (800) 945 6154 Tel: +1 512 837 7190 Fax: +1 512 836 9159 Email: usa@sge.com</p>	
<p>JAPAN SGE Japan Inc Tel: +81 45 222 2885 Fax: +81 45 222 2887 Email: japan@sge.com</p>			 <p>SGE Analytical Science www.sge.com</p>