

Nerolidol (+/- , cis - trans)
Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DET Beta - 0.25mm, 0.25µm, 25m**


Catalog Code: DEX-DET-B-025-025-25

Conditions:

Injection: Split 250°C, 0.5µL, 1:50 Split Ratio.

Sample dilution: 1mg/1mL each in n-Hexane.

Contact us for additional information.

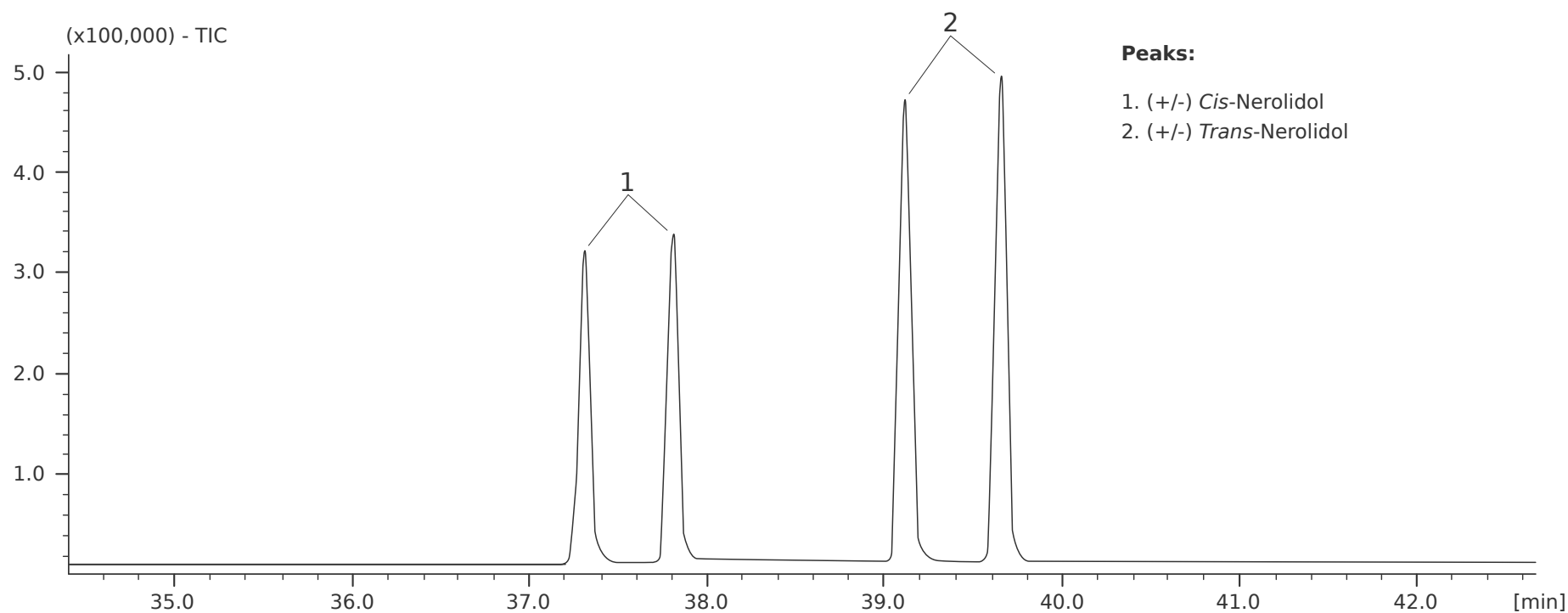

Peaks:

1. (+/-) *Cis*-Nerolidol
2. (+/-) *Trans*-Nerolidol

Acknowledgement: Prof. C.Bicchi, C.Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.

